




INTRODUCING Sonic Skim

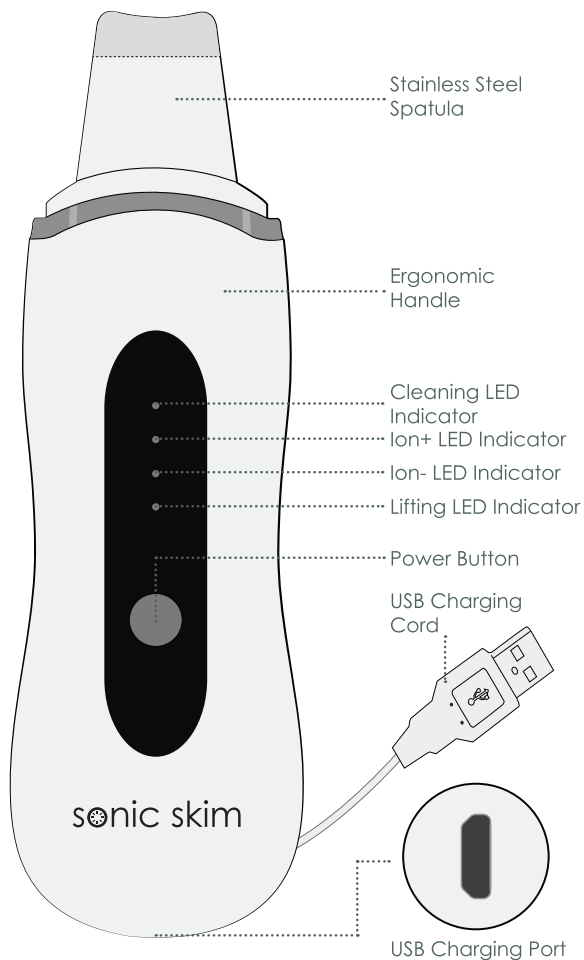


CHARGING

BEFORE FIRST USE

Sonic Skim is cordless and rechargeable. Before use, charge fully for 12 hours, even if the indicator light has changed color. Sonic Skim can run up to 3 hours on a single charge.

CHARGING

To recharge, insert the included USB Cord into the port on the side of the Sonic Skim. Then connect the USB end of the charging cord to a USB port of a computer or a USB wall adapter and then plug into an electrical outlet. As the Sonic Skim charges the white lights will illuminate from bottom to top, and once it is fully charged all the white lights will be solid white.

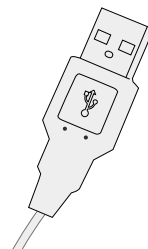

Please be sure to allow the Sonic Skim to charge for the full 12 hours, even if the white lights are solid before the 12 hours, allow it to charge for the full 12 hours. Recharge the Sonic Skim only as needed and up to 12 hours only. DO NOT overcharge which can shorten the battery life. Disconnect the Charging Cord when not charging the device. Sonic Skim should be powered OFF while charging. Do not charge Sonic Skim for more than 24 hours.

AUTOMATIC SHUT OFF

Please keep in mind that the device will shut off automatically after 15 minutes of continuous use.

HIGH FREQUENCY

Please note that the High Frequency of the Sonic Skim may be uncomfortable to pets' ears. Keep pets away while using Sonic Skim.

STEP BY STEP

READY TO USE:

Once you have assured of the full charge of your Sonic Skim device you can proceed with your treatment.

CLEANING & ION+:

It is recommended that you wash your face with cleanser and warm water before proceeding with Sonic Skim. Pull your hair back and away from your face as well. Remove the cap and press and hold the Power button for 2 seconds to power the device on - all 4 white lights will be illuminated. Press the power button one more time and this will start at the "Cleaning" Modality, press it once more when completed and you will perform the Ion+ treatment. Ion+ treatment removes the dirt and debris from your skin at a particle level after the Cleaning treatment has been complete. Make sure that you always keep your face moist with warm water when performing this treatment. Position spatula with the tip angled down at a 45-degree angle, and hold skin taut with finger from the opposite hand. Glide device gently over face and neck, while applying gentle pressure, to extract debris from pores. Visible particles/debris from your pores will build up on the spatula. Wipe clean as needed throughout your treatment.



EXPERT ADVICE

EXPERT ADVICE

WARNINGS & WARRANTY

REGISTRATION REQUIREMENT.

YOU MUST REGISTER YOUR PRODUCT WITHIN 90 DAYS OF PURCHASE IN ORDER TO OBTAIN THE WARRANTY. GO ONLINE TO WWW.WARRANTYMICHAELTODDBEAUTY.COM AND COMPLETE THE INFORMATION REQUIRED FOR REGISTRATION.

WARRANTY LIMITATIONS.

TO THE EXTENT PERMITTED BY LAW, THIS WARRANTY AND THE REMEDIES SET FORTH ARE EXCLUSIVE AND IN PLACE OF ALL OTHER WARRANTIES, REMEDIES AND CONDITIONS, WHETHER ORAL OR WRITTEN, STATUTORY, EXPRESS OR IMPLIED. MICHAEL TODD BEAUTY ("MICHAEL TODD") DISCLAIMS ALL STATUTORY AND IMPLIED WARRANTIES, INCLUDING WITHOUT LIMITATION WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE AND WARRANTIES AGAINST HIDDEN OR LATENT DEFECTS TO THE EXTENT PERMITTED BY LAW. IN ALL EVENTS, MICHAEL TODD LIMITS THE DURATION AND REMEDIES OF SUCH WARRANTIES TO THE DURATION OF THIS EXPRESS WARRANTY AND AT MICHAEL TODD'S OPTION THE REPAIR OR REPLACEMENT SERVICES DESCRIBED BELOW.

LIMITATION OF LIABILITY.

EXCEPT AS PROVIDED IN THIS WARRANTY AND TO THE MAXIMUM EXTENT PERMITTED BY LAW, MICHAEL TODD IS NOT RESPONSIBLE FOR DIRECT, SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES RESULTING FROM ANY BREACH OF WARRANTY OR CONDITION, OR UNDER ANY OTHER LEGAL THEORY, INCLUDING BUT NOT LIMITED TO LOSS OF USE; LOSS OF REVENUE; LOSS OF BUSINESS; LOSS OF OPPORTUNITY; LOSS OF GOODWILL; LOSS OF REPUTATION; OR ANY INDIRECT OR CONSEQUENTIAL LOSS OR DAMAGE HOWSOEVER CAUSED INCLUDING THE REPLACEMENT OF EQUIPMENT AND PROPERTY. THE FOREGOING LIMITATION SHALL NOT APPLY TO DEATH OR PERSONAL INJURY CLAIMS, OR ANY STATUTORY LIABILITY FOR INTENTIONAL AND GROSS NEGLIGENT ACTS AND/OR OMISSIONS.

WHAT IS COVERED BY THIS WARRANTY.

MICHAEL TODD warrants this product against defects in materials and workmanship when used normally in accordance with MICHAEL TODD'S guidelines (including technical specifications, user manuals and service communications) for a period of THREE (3) YEARS from the date of original retail purchase ("Warranty Period"). This warranty applies only to a purchase of the Product in the United States or Canada from MICHAEL TODD or an authorized reseller, is valid for the original end user of the Product only and cannot be transferred under any circumstances.

NEED MORE INFO?

The entire operating manual including safety, care and warranty information is available at

www.sonicskimmanual.com

Contact customer service at 772-800-7654 or customerservice@michaeltoddbeauty.com

© 2021 MICHAEL TODD BEAUTY LP. All rights reserved.

QUICK TIPS TO OPTIMIZE YOUR RESULTS WITH SONIC SKIM:

1. Sonic Skim is equipped with Sonic Technology and Medical-Grade Stainless Steel treatment Spatula for Face and Body Treatments and designed to replicate the spiral kneading action of an Esthetician's fingers.
2. The sonic technology has more than 30,000 movements per minute. It gently massages, tones your skin and increases circulation.
3. Treatment Spatula for Face and Body Treatments and designed to replicate the spiral kneading action of an Esthetician's fingers.
4. This device is negatively charged device that helps to push serums/creams deeper into the skin with the ION-technology.
5. Treatment Spatula for Face and Body Treatments and designed to replicate the spiral kneading action of an Esthetician's fingers.
6. This device negatively charged function helps to push serums/creams deeper into the skin with the ION-technology.
7. Use the Lifting modality, alone or in combination with the other functions of the Sonic Skim. The EMS (Electro-Muscle Stimulation) technology helps tighten and lift your skin for a tighter more youthful appearance.

IMPORTANT:

8. DO NOT dig the corners of the Sonic Skim into your skin as this can damage your skin causing redness. Sonic Skim is to be used on uncompromised treatment areas. Avoid use in treatment areas where any of the following is present; any area that was recently treated with chemical peels or lasers, over raised areas of skin such as lesions, moles, acne, and skin tags. Avoid treatment areas where cancer is present or suspected to be present. DO NOT use over compromised skin. Compromised skin may include but is not limited to sunburns, areas of healing, blisters, Rosacea, Psoriasis, abrasions, cuts and acne.

BASIC SAFETY PRECAUTIONS.

When using electrical products, especially when children are present, basic safety precautions should always be taken. IF YOU THINK THE POWER IS OFF WHEN THE SWITCH IS OFF. YOU ARE WRONG. If the plug is in, the power is on. Unplug electrical appliances when you are not using them.

DANGERS.

To reduce the risk of electrocution: DO NOT get the product or Charging Cable wet or immerse in water. Keep this product away from open flame or flammable aerosol products. The Device is not suitable for use in the presence of a flammable anesthetic mixture with air or oxygen or nitrous oxide.

WARNINGS.

To reduce the risk of burns, electrocution, fire or injury to persons: Close supervision is necessary when this product is used by, on, or near children or physically challenged persons. Use this product only for its intended use as described in this User Guide. DO NOT use this product where oxygen is being administered. The product has no serviceable parts. Never attempt to open the Device Handle for any reason. DO NOT disassemble, heat above 600°C (140o°F) or incinerate.

IMPORTANT CAUTIONS.

This is not a Medical Device and is not intended to treat medical conditions. If you are under a doctor's care for any skin condition. if you think you may have a skin condition or if you have very sensitive skin, consult your physician before using this product. This product uses magnets which may affect internal medical devices such as pacemakers. Keep away from internal medical devices. This Device is electrical and contains small parts. Keep out of reach of children. THIS PRODUCT IS NOT DESIGNED OR INTENDED PRIMARILY FOR CHILDREN 12 YEARS OF AGE OR YOUNGER.

If you have any medical concerns about using Sonic Skim, please consult your dermatologist before use.

If any of the components are damaged (Sonic Skim Device Handle, or USB Charging Cord) DO NOT use. DO NOT attempt to open or service components. Email Michael Todd Beauty customer services at customerservice@michaeltoddbeauty.com

ION-:

Remove the cap and press and hold the Power button for 2 seconds to power the device on all 4 white lights will be illuminated. Press the power button three more times and this will start at the "ION-" Modality. Apply your serums/cream or moisturizers to your skin. Use this with the tip facing out away from your skin, using the opposite side you used from the Cleaning mode. The negatively charged Sonic Skim will help to push your products about 3-5 mm into the skin for maximum benefits of your products. Drag your Sonic Skim slowly and making multiple passes, which will help to maximize the product infusion into your skin.

LIFTING:

To tighten and lift your skin, press and hold the power button to activate the device and press four times on the power button until the fourth white light is illuminated. Place the device on your skin with the tip facing up and away from your skin, using the same side you used in the Ion - mode. Drag your Sonic Skim slowly and making multiple passes over the same area to help enhance the lifting and tighten effect that the EMS provides to your skin.

CLEANING OF THE SPATULA:

You may wipe your Sonic Skim's Stainless Steel Spatula at the end of your treatment with a microfiber cloth. Use exclusively non-abrasive cleaning liquids such as warm water, mild dish soap, or a solution with part of distilled white vinegar or alcohol. Always spot test before using on the whole Spatula. NEVER submerge the whole device/handle in any liquid. NEVER use steel wool, bleach, scouring powders or ammonia. In the long run, hard water can also leave spots or streaks.